

FILTER BARRIER TIPS

Molded in virgin medical grade polypropylene, NOVASBIO filter barrier tips provide exceptional performance when pipetting critical solutions.




Aerosol barrier tips - prevent aerosol and liquid contamination

Available in Low Retention for complete transfer of valuable samples

Available in bulk and racked formats

Pre-sterilized by E-beam irradiation

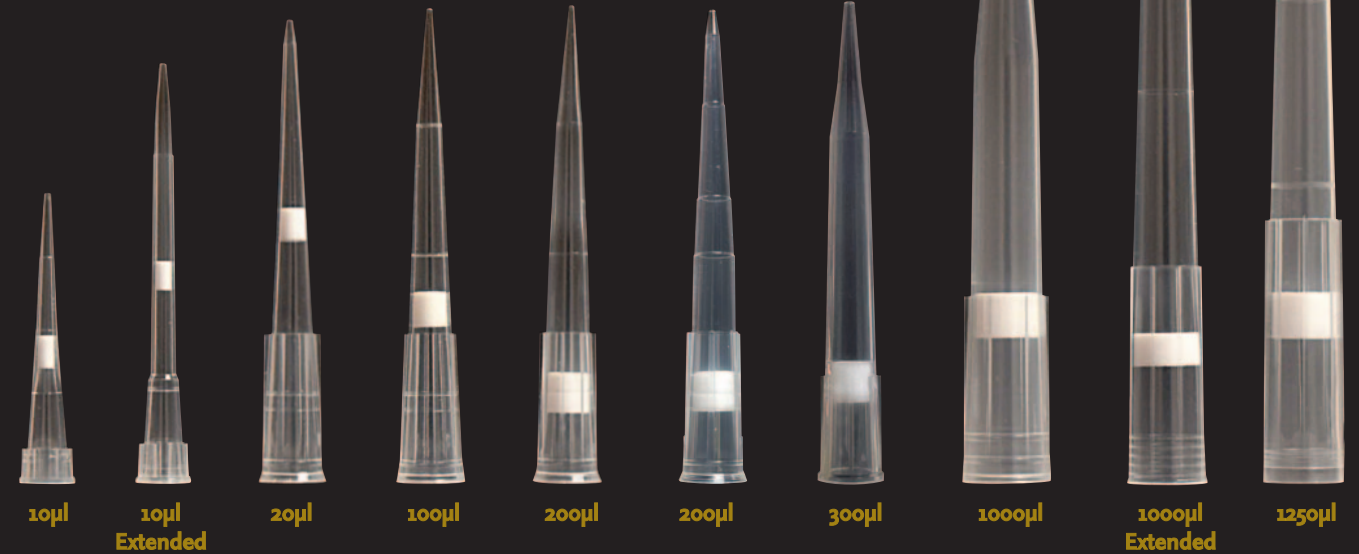

The thin-wall design offers a perfect seal and excellent compatibility with most popular brands of pipettors.

The rated volume of each tip can be aspirated with space to spare, while the inert Ultra High Molecular Weight Polyethylene (UHMWPE) filter with its consistent porosity (20 microns) allows air flow but not aqueous solutions, so protecting samples and pipettor from contamination.



RACKED TIPS
Clasp-lock keeps the lid securely fastened while allowing one-handed operation.

Low Retention

NOVASBIO's low retention manufacturing process results in an ultra-hydrophobic surface, preventing liquid adhering to the surface and so ensuring maximum sample transfer. This NOVASBIO low retention process is implemented pre-production and does not involve surface coatings that could lead to sample contamination.

CATALOG NO.	VOLUME	DESCRIPTION	TIP LENGTH	PACK	CASE
FT-10-R-S	10µl	Filtered, Racked, Graduations, Pre-Sterilized	31mm	10 x Rack/96 Tips	5 Packs
FT-10-L-R-S	10µl	Filtered, Low Retention, Racked, Graduations, Pre-Sterilized	31mm	10 x Rack/96 Tips	5 Packs
FT-10XL-R-S	10µl	Extended, Racked, Pre-Sterilized	45mm	10 x Rack/96 Tips	5 Packs
FT-10XL-L-R-S	10µl	Extended, Low Retention, Racked, Pre-Sterilized	45mm	10 x Rack/96 Tips	5 Packs
FT-20-R-S	20µl	Filtered, Racked, Pre-Sterilized	51mm	10 x Rack/96 Tips	5 Packs
FT-20-L-R-S	20µl	Filtered, Low Retention, Racked, Pre-Sterilized	51mm	10 x Rack/96 Tips	5 Packs
FT-100-R-S	100µl	Filtered, Racked, Pre-Sterilized	52mm	10 x Rack/96 Tips	5 Packs
FT-100-L-R-S	100µl	Filtered, Low Retention, Racked, Pre-Sterilized	52mm	10 x Rack/96 Tips	5 Packs
FT-200-R-S	200µl	Filtered, Racked, Pre-Sterilized	52mm	10 x Rack/96 Tips	5 Packs
FT-200-L-R-S	200µl	Filtered, Low Retention, Racked, Pre-Sterilized	52mm	10 x Rack/96 Tips	5 Packs
FT-201-L-R-S	200µl	Filtered, Low Retention, Racked, Graduations, Pre-Sterilized	51mm	10 x Rack/96 Tips	5 Packs
FT-300-L-R-S	300µl	Filtered, Low Retention, Racked, Pre-Sterilized	52mm	10 x Rack/96 Tips	5 Packs
FT-1000-L-R-S	1000µl	Filtered, Low Retention, Racked, Pre-Sterilized	78mm	10 x Rack/100 Tips	5 Packs
FT-1000XL-R-S	1000µl	Filtered, Extended, Racked, Graduations, Pre-Sterilized	83mm	10 x Rack/96 Tips	5 Packs
FT-1000XL-L-R-S	1000µl	Filtered, Ext., Low Retention, Racked, Graduations, Pre-Sterilized	83mm	10 x Rack/96 Tips	5 Packs
FT-1250-L-R-S	1250µl	Filtered, Ext., Low Retention, Racked, Graduations, Pre-Sterilized	102mm	8 x Rack/96 Tips	5 Packs

* on 10µl, 100µl, 200µl, 1000µl Extended and 1250µl tips only. Bulk, non-sterile filter tips available upon request.

