

NOVASBIO Screw Cap Microtubes are used in a wide range of laboratory applications. These include collection, processing, centrifugation, long-term storage and packaging/shipping of valuable biological materials.

Ideal for long-term, secure storage and transportation of biological reagents

Available with 5 colors of caps plus clear


Suitable for high speed centrifugation up to 20,000 x g

Can be used at extreme temperatures ranging from -80°C to +121°C \*



Color-coded caps provide a syringe port which is easily pierced for access to the contents. The threaded 'O' ring design ensures a uniform closure allowing centrifugation up to 20,000 RCF. Ribs are designed to facilitate single-handed opening and closing, even with gloves on.


NOVASBIO screw cap microtubes are supplied non-sterile, un-assembled with both Tubes and Caps packed separately. To special order, they are available pre-sterilized, assembled and packed in convenient bags of 100.

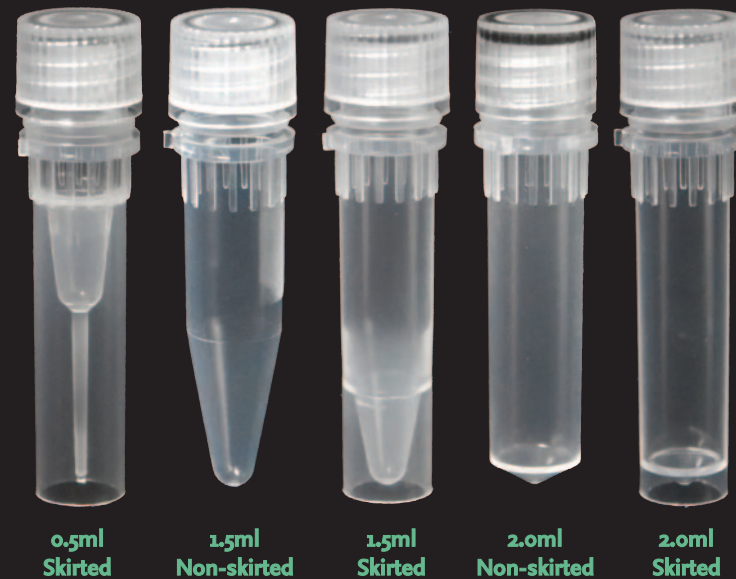


CATALOG NO.	VOLUME	DESCRIPTION	0.5ml
SCT-050SS-C	0.5ml	Self Standing Clear Tube, Clear Screw Cap, Non Sterile	
SCT-050SS-G	0.5ml	Self Standing Clear Tube, Green Screw Cap, Non Sterile	
SCT-050SS-W	0.5ml	Self Standing Clear Tube, White Screw Cap, Non Sterile	
SCT-050SS-B	0.5ml	Self Standing Clear Tube, Blue Screw Cap, Non Sterile	
SCT-050SS-Y	0.5ml	Self Standing Clear Tube, Yellow Screw Cap, Non Sterile	
SCT-050SS-R	0.5ml	Self Standing Clear Tube, Red Screw Cap, Non Sterile	




CATALOG NO.	VOLUME	DESCRIPTION	1.5ml
SCT-150-C	1.5ml	Conical Bottom Clear Tube, Clear Screw Cap, Non Sterile	
SCT-150-G	1.5ml	Conical Bottom Clear Tube, Green Screw Cap, Non Sterile	
SCT-150-W	1.5ml	Conical Bottom Clear Tube, White Screw Cap, Non Sterile	
SCT-150-B	1.5ml	Conical Bottom Clear Tube, Blue Screw Cap, Non Sterile	
SCT-150-Y	1.5ml	Conical Bottom Clear Tube, Yellow Screw Cap, Non Sterile	
SCT-150-R	1.5ml	Conical Bottom Clear Tube, Red Screw Cap, Non Sterile	

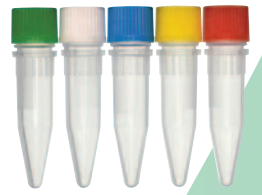


SCT-150SS-C	1.5ml	Self Standing Clear Tube, Clear Screw Cap, Non Sterile	1.5ml
SCT-150SS-G	1.5ml	Self Standing Clear Tube, Green Screw Cap, Non Sterile	
SCT-150SS-W	1.5ml	Self Standing Clear Tube, White Screw Cap, Non Sterile	
SCT-150SS-B	1.5ml	Self Standing Clear Tube, Blue Screw Cap, Non Sterile	
SCT-150SS-Y	1.5ml	Self Standing Clear Tube, Yellow Screw Cap, Non Sterile	
SCT-150SS-R	1.5ml	Self Standing Clear Tube, Red Screw Cap, Non Sterile	

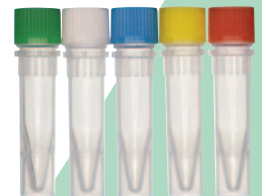

CATALOG NO.	VOLUME	DESCRIPTION	2.0ml
SCT-200-C	2.0ml	Conical Bottom Clear Tube, Clear Screw Cap, Non Sterile	
SCT-200-G	2.0ml	Conical Bottom Clear Tube, Green Screw Cap, Non Sterile	
SCT-200-W	2.0ml	Conical Bottom Clear Tube, White Screw Cap, Non Sterile	
SCT-200-B	2.0ml	Conical Bottom Clear Tube, Blue Screw Cap, Non Sterile	
SCT-200-Y	2.0ml	Conical Bottom Clear Tube, Yellow Screw Cap, Non Sterile	
SCT-200-R	2.0ml	Conical Bottom Clear Tube, Red Screw Cap, Non Sterile	

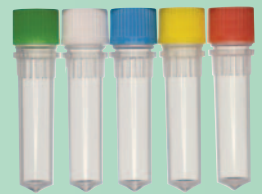

SCT-200SS-C	2.0ml	Self Standing Clear Tube, Clear Screw Cap, Non Sterile	2.0ml
SCT-200SS-G	2.0ml	Self Standing Clear Tube, Green Screw Cap, Non Sterile	
SCT-200SS-W	2.0ml	Self Standing Clear Tube, White Screw Cap, Non Sterile	
SCT-200SS-B	2.0ml	Self Standing Clear Tube, Blue Screw Cap, Non Sterile	
SCT-200SS-Y	2.0ml	Self Standing Clear Tube, Yellow Screw Cap, Non Sterile	
SCT-200SS-R	2.0ml	Self Standing Clear Tube, Red Screw Cap, Non Sterile	

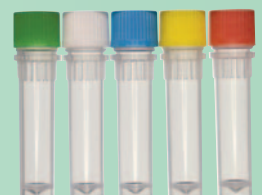


Screw Cap Tubes are also available pre-sterilized to special order, 100 Tubes per Unit, 5 Units per Pack, 2 Units per Case, Assembled. Please call for more information

**Packaging Description:** 500 per Pack, 2 Pack per Case, Unassembled

**Maximum G-force:** 20,000 RCF



**Light blocking Amber**

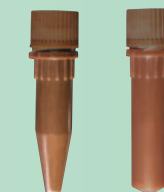

Available in light blocking Amber with Brown Screw Caps upon special order